

Instructions for Attaching a Stand Accessory to a 1 Series Phototherapy Device

Tools Needed: A Phillips-head Screwdriver (has “+” shaped end)

1 Series Stand Accessory Kit Includes: 1 Cradle, 2 Spacer Disks, 2 Knobs, and 4 Phillips Screws

Step 1: Unplug the power cord from the device.

Step 2: Find a suitable work space and stand the device upright so that the lamps are facing toward you and the timer is towards the top, shown in the illustration at right.




Step 3: Place one of the round spacer disks on the top of the device and align the two holes on the disk with the two holes in the device exterior.

Step 4: Insert two screws into the holes and tighten with a Phillips screwdriver.

Step 5: Place the other round spacer disk on the bottom side of the device and align with the holes there.

Step 6: Insert two screws into the holes and tighten with a Phillips screwdriver.

Step 7: Place the cradle next to the device and align the holes on the cradle with the holes on the spacer disks.

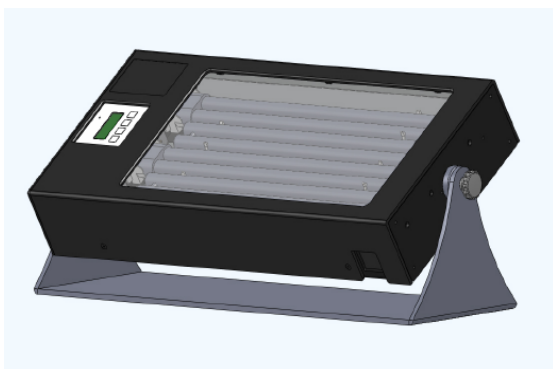

Step 8: Insert the black knobs into the holes and screw in to tighten. Your device is now ready for use, as shown in the illustration at left.